The server core (**dgate.exe** = **dgate** under Linux) compiles and runs on Linux systems and Solaris. I develop primarily under Windows, but currently I test the code and scripts under Linux Ubuntu 18.04. I also had the server compiled on a Raspberry Pi.

The Linux release of the server core works default with SqLite driver built in into the server (no ODBC). The DbaseIII driver is also supported. Piotr Filipczuk has added a PostGresQL driver. The native MySQL interface also can be used. The graphical user interface has not been ported to Linux, but the WEB interface is provided, either using Apache or a built-in mini web server (Ladle). In this version, most options have been well tested – it is a stable release. However, there are often subtle differences between linux distributions, making installation (and writing a manual) difficult. There are several contributions on the forum, and there are text files with specific command orders to be found in the linux subfolder of the server.

To use the server, one needs a valid version of the configuration files and put them in the same directory as the dgate executable. The easiest way to do this is to unpack **dicomserver150beta3.zip** with "unzip dicomserver150beta3.zip".

INSTALLATION

Prerequisites: 1) a running Linux system. 2) sudo installed and enough rights. If not, the script will not be able to install the server as web service for Apache and you need to copy the files by hand. Note that I only test the scripts on Ubuntu, but the web based installer script linux.sh has a bit of info on Fedora.

These packages needed to be installed in a plain Linux system for a release using SQLite or DbaseIII:

```
sudo apt update
                                                                                          get compilers
sudo apt install make
                                                                                          may be skipped if using precompiled
(or: sudo apt install build-essential)
sudo apt install q++
                                                                                          may be skipped if using precompiled
sudo apt install apache2
                                                                                          aet webserver
sudo apt install unzip
                                                                                          not standard in Ubuntu server
                                                                                          parts of the web interface use 7za
sudo apt install p7zip-full
sudo apt install lua5.1
                                                                                          since 1.5.0beta3 lua is external
                                                                                          may be skipped if using precompiled
sudo apt install lua5.1-dev
(or sudo apt install liblua5.1-0)
                                                                                          is required when using precompiled
sudo apt install lua-socket
(or for fedora:
dnf install gcc-c++-sh-linux-gnu.x86 64 gcc-c++-x86 64-linux-gnu.x86 64 clang.x86 64
sudo a2enmod cgi
                                                                                          enable CGI in web server
systemctl restart apache2
(or for older systems: sudo service apache2 restart)
```

The rest of the installation can be performed manually, or by a web based method, explained below.

The following steps illustrate a minimal installation:

First get the server:

wget http://ingenium.home.xs4all.nl/dicomserver/dicomserver150beta3.zip

get server zip

mkdir conquest cd conquest

unzip ../dicomserver150beta3.zip rm ../dicomserver150beta3.zip

make folder to store conquest to there

Or:

sudo apt install git git clone https://github.com/marcelvanherk/Conquest-DICOM-Server cd Conquest-DICOM-Server

if git not installed yet get latest from GitHub

Then compile and install it:

chmod 777 maklinux ./maklinux choose option 3 or 5 say 'y' to 'Regenerate the database' say 'y' to 'Install as service' Shows status hit 'q' to return

compile and install web access Salite or Sallite precompiled

Deletes previous database contents

Now the server should be running and http://localhost/cgi-bin/newweb/dgate should provide a working web interface

Note that in dicomserver150beta3 a precompiled dgate (compiled by me on Ubuntu 14.04, using Sqlite database) is included, to try that use option 5 in maklinux. Tested on Ubuntu 18.04, 19.10. If used the following packages may be *omitted*: **make**, **g++**, **lua5.1-dev**; but if you do omot them then the following package must be added: **liblua 5.1-0.** This option reduces the size of the Linux system by a few hundred MB.

Web based installation

To run the web based installer (after installing prerequisites):

wget http://ingenium.home.xs4all.nl/dicomserver/dicomserver150betae.zip

mkdir conquest

cd conquest

unzip ../dicomserver150beta3.zip

cd install

chmod 777 linux.sh

./linux.sh


get server zip

make folder to store conquest

make run-able

This compiles a minimal server binary (dgatesmall) or uses the pre-compiled one, that is run as service

control manager and, if a web server and client (Firefox expected) exist, opens web page http://127.0.0.1/cgi-bin/service/dgate. The resulting web page allows and guides the user through compilation, configuration, re-generation of the database if needes, starting the server, setting up the web server and opening the web client. A screen-shot of the install page is shown below:

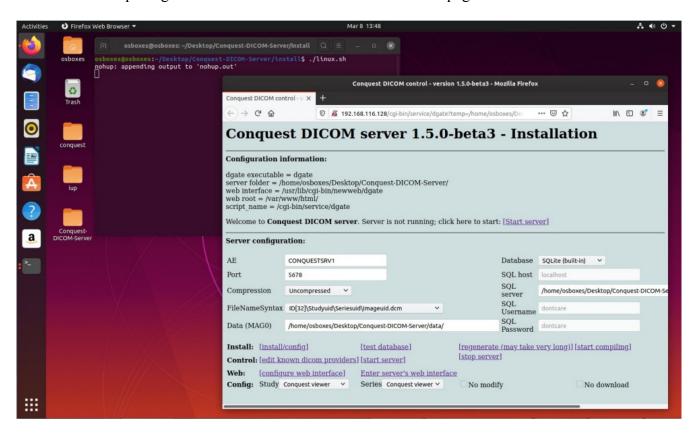

The required steps (most are shown in the welcome area) are:

- 1) Select required database type (start with SQLite if unsure)
- 2) Start compiling → compile jpeg6c, compile openjpeg, compile charls, compile lua, compile luasocket, compile sqlite, compile dgate; [done].
 - If any of the compilation steps fails error messages can be found in file nohup.out. If the compilation information disappears click start compiling again.
- 3) Set other parameters (keep defaults if unsure)
- 4) Configure server
- 5) Start server (may need be repeated a few times if does not start)
- 6) Regenerate database
- 7) Configure web interface (select viewers and access rights)

 Note that if you want to use weasis as viewer, the weasis folder from weasis_portable.zip (v3), must be placed in your web servers root folder
- 8) Enter server's web interface

Feedback on this new installation method would be appreciated. After installation, the server runs as part of the control manager. To make it run permanently, stop the server control manager (dgatesmall) with ^C, and use the new start-stop-daemon method described above or the old one below. Note that stopping the server using this web page on Linux disables restarting it for a minute or so (due to an IP port being blocked). Be patient when it fails not restart and try again after a while.

Deamon configuration

Both the web install and malinix now create a daemon as follow, changing the file to point to the conquest installation:

sudo cp conquest.service /etc/systemd/system/conquest.service sudo systemctl daemon-reload

After installation you can control the conquest service as follows:

sudo systemctl start conquest.service sudo systemctl enable conquest.service sudo systemctl status conquest.service

hit 'Q' to return

sudo systemctl stop conquest.service sudo systemctl disable conquest.service


Now the server should be running, also after a system restart, and localhost/cgi-bin/newweb/dgate should provide a working web interface.

Built-in WEB Viewer

A new single user web viewer can be run as follows:

chmod 777 webviewer.sh /webviewer.sh

This is the same web viewer as can be accessed from a full featured web server, but instead it runs on 127.0.0.1:8086, using Ladle (single user web server) as mini web server. After stopping the browser, the Ladle function is stopped. It takes a minute or so for the used port (8086) to be released. Until then attempting to start the web viewer fails.

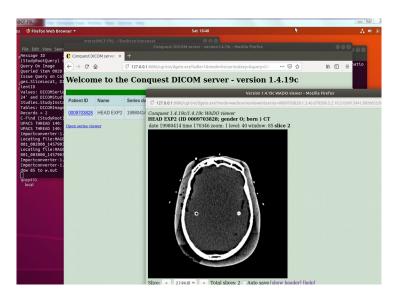

Installing with Postgres

To install with Postgres as database, these commands are needed to install and setup Postgres:

sudo apt-get install libpg-dev sudo apt-get install postgresgl sudo su

su - postgres lpzq

\password

postgres postgres

\q

createdb conquest

exit exit

./maklinux choose option 2

Series Database

Postgres development tools

Postgres database become superuser

become postgres user

set the passwork to postgres

(password) (repeat password)

create database conquest

compile and install web access

Postgres

The build process always gives a few error messages that can be ignored:

/usr/bin/install: cannot create regular file '/usr/local/man/man1/cjpeg.1': No such file or directory Makefile:200: recipe for target 'install' failed mkdir: cannot create directory 'data/dbase': File exists

During database creation (dgate -v -r) there can be error messages about non-existing databases, e.g. for postgress:

```
osboxes@osboxes:~/Desktop/distribution$ ./dgate -v -r
Regen Database
Step 1: Re-intialize SQL Tables
*** ERROR: relation "dicomworklist" does not exist
***Error: ERROR: table "uidmods" does not exist
WorkList Database
Patient Database
Study Database
```

Image Database


Step 2: Load / Add DICOM Object files

Regen Device 'MAGO'

[Regen] ./data/0009703828/1.3.46.670589.5.2.10.2156913941.892665339.860724 0001 003000 14579035620001.dcm -SUCCESS

Regeneration Complete

osboxes@osboxes:~/Desktop/distribution\$./dgate -v

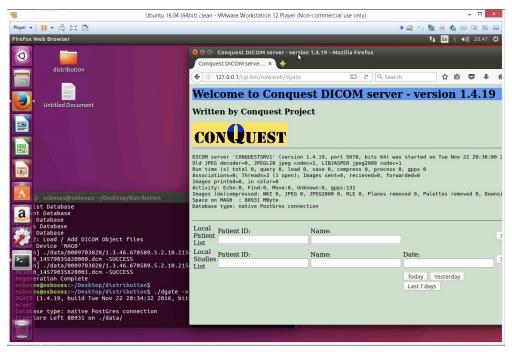

Conquest in action on Ubuntu16.04, with Postgres database and web interface

Installing with Mariadb

To install with Mariadb as database, these commands are needed to install and setup:


Mariadb server sudo apt install mariadb-server sudo apt install libmariadbolient-dev Client code sudo mysgl database superuser >create user conquest; >grant all privileges on *.* to conquest@localhost identified by 'conquest'; create user >create database conquest; create database >flush privileges; **>**\q ./maklinux compile and install web access choose option 1 mariadb

ZerobraneStudio IDE

To install and use ZeroBrane Studio with the conquest DICOM server under Linux, take these steps. First download ZeroBraneStudioEduPack-xxx-linux.sh. Then in a command prompt run:

chmod 777 ZeroBraneStudioEduPack-xxx-linux.sh sudo ./ZeroBraneStudioEduPack-xxx-linux.sh

After installation is done run ZeroBrane Studio from the command prompt as "sudo zbstudio" and run the install script /dicomserver/ZeroBraneStudio/install.lua in ZeroBrane Studio as described in this file. After running the conquest install script as root, ZeroBraneStudio can be run as a normal user.

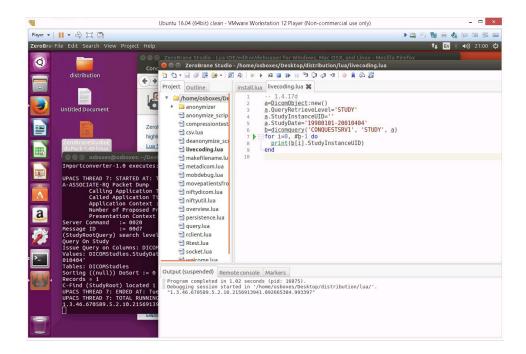

Integration with Zerobrane studio

CONFIGURATION

Configuration files under Windows and Linux are the same except for the use of a forward slash instead of back slash in directory paths. The following essential entries are therefore different for Linux (these are the defaults):

SQLServer = ./data/dbase/conquest.db3

MAGDevice0 = ./data/

See the Windows manual for more details about the configuration files (you need at least to edit acrnema.map to define DICOM systems that will be retrieving information from your server). All configurations options in dicom.ini (e.g., for DICOM routing) are listed in windowsmanual.pdf. You probably also need to edit the web server configuration file /usr/lib/cgi-bin/dicom.ini to set the correct IP address of the machine. If not the web server will only partly function.

After copying the files, if needed, regenerate the database with "conquest/dgate -v -r" then run the server with "conquest/dgate -v B" or "conquest/dgate -^serverstatus.log". NOTE: regeneration is only needed after an upgrade if **dicom.sql** is updated. If you want to avoid regeneration do NOT replace **dicom.sql**

The building process for the server was tested with gcc 3.3.5, Ubuntu 8.10 and on Solaris 10. Both 32 and 64 bit OS's are supported. Warnings are produced but these do not impact server operation.

Also MySQL support is provided. It requires creating a DB called "conquest" with phpmyadmin and installing libmysqlclientdev with: "apt-get install libmysqlclient-dev" before running maklinux_mysql. These are the settings in dicom.ini for MySQL:

SQLHost = localhost
SQLServer = conquest
Username = root
Password =
Mysql = 1
DoubleBackSlashToDB = 1

The PostGres system can be setup to the defaults, and a database named 'conquest' made. For postgres to work you need to check some values in dicom.ini (using the default postgres account assuming password postgres, note that parameter 'SQLServer' sets the database to conquest). A copy from dicom.ini.postgres to dicom.ini would set the following values:

SQLHost = localhost SQLServer = conquest Username = postgres Password = postgres

PostGres = 1
DoubleBackSlashToDB = 1
UseEscapeStringConstants = 1

It is advised to use a normalized database (as defined in **dicom.sql**) for postgres operation, e.g., by copying **dicom.sql.postgres** to **dicom.sql** and a denormalized database for DbaseIII, e.g., by copying **dicom.sql.dbase** to **dicom.sql**.

The following are donated scripts by Mark Pearson for start/stop and rotating logfiles and this information is for expert users only:

To install this script (it is in the distribution as nconquest-pacs.sh) do:

sudo cp nconquest-pacs.sh /etc/init.d/ sudo chmod 755 /etc/init.d/nconquest-pacs.sh sudo apt-get install authbind sudo /etc/init.d/nconquest-pacs.sh start

```
#!/bin/bash
                        SysV init script for Conquest PACS.
  conquest-pacs.sh
        Written by Miquel van Smoorenburg <miquels>.
        Modified for Debian GNU/Linux by Ian Murdock <imurdock>.
        Customized for Conquest by Mark Pearson <markp>
        HOME and PACSUSER should be the only variables that may need to be modified.
PATH=/sbin:/bin:/usr/sbin:/usr/bin
# Modify HOME to suit your environment.
HOME=/usr/local/conquest
# This is the user to run as. Modify it if you don't use username conquest.
PACSUSER=conquest
DAEMON=$HOME/dgate
INI=$HOME/dicom.ini
NAME=conquest_pacs.sh
# All defaults here will be overridden by values from $HOME/dicom.ini
STATUSLOG=$HOME/serverstatus.log
DESC="Conquest PACS Server"
STOPPACS=$HOME"/dgate --quit:"
STARTAS=$DAEMON
test -f $DAEMON || echo "Cannot find $DAEMON" exit 0
test -f $INI || echo "Cannot find $INI" exit 0
set -e
if grep "TCPPort" $INI > /dev/null; then
        PORT=`egrep -i '^*TCPPort *= ' $INI | sed 's/\r//' | awk '{ print $3}'`
fi
if [ $PORT -le 1024 ]; then
        test -f /usr/bin/authbind || echo "authbind is needed for access to ports < 1024" exit 0
        STARTAS="/usr/bin/authbind "
fi
if grep -is "^ *StatusLog" $INI > /dev/null; then
        STATUSLOG=`egrep -i '^*StatusLog' $INI | sed 's/\r//' | awk '{ print $3}'`
fi
PIDFILE=/var/run/$NAME.$PORT.pid
if [ $STARTAS = $DAEMON ]; then
       ARGS=" -^$STATUSLOG"
        ARGS="$DAEMON -^$STATUSLOG"
fi
case "$1" in
        if [ -f $HOME/disable_autostart ]; then
                echo "Not starting $DESC: disabled via $HOME/disable autostart"
                exit 0
        fi
```

```
echo -n "Starting $DESC: "
        start-stop-daemon --start --quiet --pidfile $PIDFILE \
                --chuid $PACSUSER --chdir $HOME --exec $DAEMON \
                --startas $STARTAS --background -- $ARGS
        echo "$NAME."
        ; ;
 stop)
        echo -n "Stopping $DESC: "
        cd $HOME
        $STOPPACS
        start-stop-daemon --oknodo --stop --quiet --pidfile $PIDFILE \
                --exec $DAEMON -- $ARGS
        echo "$NAME."
        echo
        ;;
  restart|force-reload)
        echo -n "Restarting $DESC: "
        start-stop-daemon --stop --oknodo --quiet --pidfile $PIDFILE \
                --exec $DAEMON -- $ARGS
        start-stop-daemon --start --quiet --pidfile $PIDFILE \
                --chuid conquest --chdir $HOME --exec $DAEMON -- $ARGS
        echo "$NAME."
  *)
        N=/etc/init.d/$NAME
        echo "Usage: $N {start|stop|restart|force-reload}" >&2
        exit 1
esac
exit 0
```

For security reasons I have added a user "conquest" and the package authbind to allow access to priveleged ports. I added the following entries to dicom ini:

HomeDir = /usr/local/conquest

 $StatusLog = \frac{1}{\sqrt{NMPACS}}. server status.log$

TroubleLog = /var/log/conquest/NMPACS.PacsTrouble.log

The file /etc/cron.weekly/conquest rotate does weekly log rotation for me.

```
#!/bin/bash
# conquest rotate
                        Cron script to rotate conquest log files.
      Keep files for 365 days
      Read filenames from dicom.ini
                 Written by Mark Pearson 20070711 <markp>.
# Modify this line to suit your environment
HOMES=(/usr/local/conquest /usr/local/conquest-icon)
for i in ${HOMES[@]}; do
        INI=${i}/dicom.ini
        STATUSLOG=${i}/serverstatus.log
        TROUBLELOG=${i}/PacsTrouble.log
        set -e
# defaults will be overridden by values from ${i}/dicom.ini
        if grep -is "^ *StatusLog" $INI > /dev/null; then
                 \label{eq:statuslog} $$\operatorname{STATUSLOG=`egrep -i '^*Statuslog' $INI \mid sed 's/\r'/' \mid awk '{ print $3}'`} $$
        fi
            grep -is "^ *TroubleLog" $INI > /dev/null ; then
        if
                 TROUBLELOG=`egrep -i '^*TroubleLog' $INI | sed 's/\r//' | awk '{ print $3}'`
        fi
```

This copes with multiple pacs instances on the same host. The advantage of using savelog is that old logfiles are compressed. It should be quite simple to edit the files to have executable or log in /opt. Also, don't forget to set the appropriate file permissions for the user that runs conquest.

Finally, Here are the command lines to compile the server under OS X xcode using 10.4u sdk on a PowerPC (not recently tested):

And to compile under SOLARIS 10:

 $/usr/sfw/bin/g++ -DUNIX -DNATIVE_ENDIAN=1 -DHAVE_LIBJASPER -DHAVE_LIBJPEG -DSOLARIS \ total.cxx -o \ dgate -lpthread -lsocket -lnsl -lposix4$